# Vitek®

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PERSONAL SCALE AND BODY FAT ANALYSER

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## ENGLISH

### BODY FAT ANALYSER

This analyser will be a help to a healthier lifestyle. Read this manual careful and keep it for future reference.

### WHAT IS BODY FAT?

It is a well known fact that increased body fat can be detrimental for your health. Too much fat can lead to heart problems, arteriosclerosis and obesity while too little fat can cause osteoporosis especially in women. However through exercise and a contained diet you can control your body fat count. Body Fat Analyser is ideal for keeping check of not only your weight but also for monitoring the percentage of fat in your body. We know that muscles and body tissues with a low body fat count contain more water and that an electrical current passes more quickly through water than through fat. By allowing a weak current to pass through the body you can measure quite simply the percentage of fat retained in your muscles. That way you have a quick and easy way of keeping check of your fitness.

### There are some things you should think about!

- When you want a body fat reading, you must have bare feet.
- Always measure yourself at the same time of day, preferably in the afternoon.
- Wait 2-3 hours after meals.
- Wait for up to 12 hours after a strenuous workout.
- Measure yourself after a visit to the bathroom.

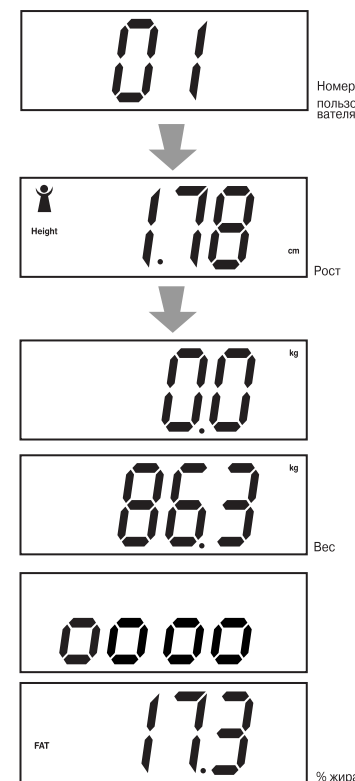
### Precautions

- **Pregnant women should NOT use a body fat analyser.**
- **Persons with pacemakers or other battery driven medical devices should NOT use the body fat analyser.**

### PROBLEMS

1. In case of abnormal use of the unit, or the body fat percentage is out of range, the display will show "Err".
2. "FFFF" will be shown when the scale is overloaded.
3. If the display shows "batt", it indicates that a new battery should be inserted.

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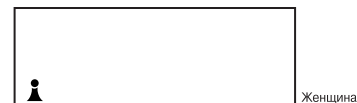
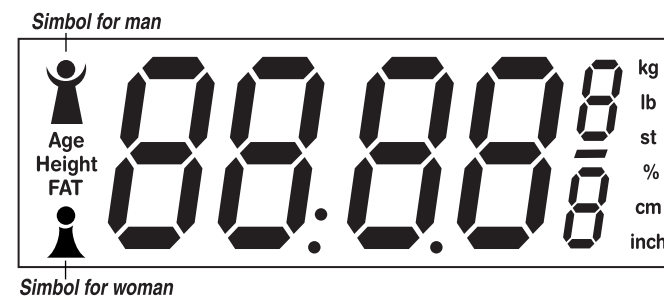
## ENGLISH

**DISPLAY****HOW TO USE THE ANALYSER****BATTERY INSTALLATION**

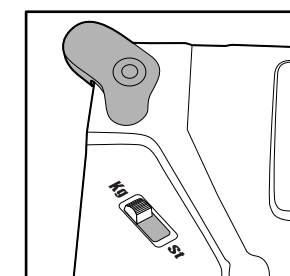
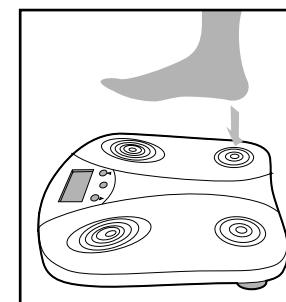
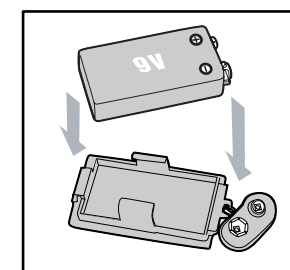
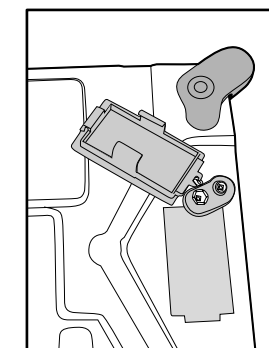
Place one new 9 volt alkaline battery (not supplied) into the battery compartment at the back of the unit.

**TO CHOOSE KILOS OR STONES**

On the back of the scale there is a small switch. If you put the switch to "KG" the measurements will be in kilos and centimeters for height. If you put the switch to "ST" the measurements will be in stones / lbs for weight and feet/inches for height.

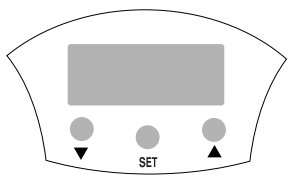
**TO TURN ON THE UNIT**

Press the bottom right hand corner of the unit using your foot.

## ENGLISH

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**USE OF KEYS**

1. SET Press to confirm data entry
2. ▲ Press to increase the number displayed
3. ▼ Press to decrease the number displayed

**WEIGHT MEASUREMENT**

1. Press the unit on with your foot.
2. The display will show "8888" and then "0.0".
3. Stand steadily on the unit until your correct weight is shown.

The weight will display for a few seconds before turning off automatically.

Weight

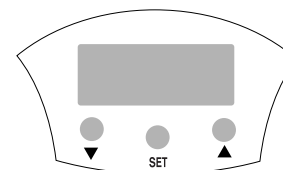
User number

**INPUT OF PERSONAL DATA**

1. Press the unit on and the display will show "8888" and then "0.0"
2. Press the SET key and the user number is displayed. Press ▲ to increase the user number or press ▼ to decrease the user number. There are four user numbers available - 01, 02, 03 and 04. Press SET to confirm.

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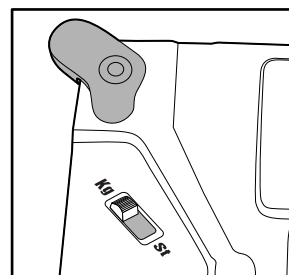
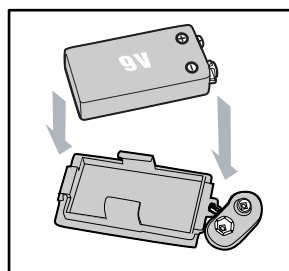
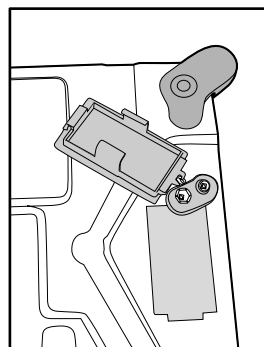
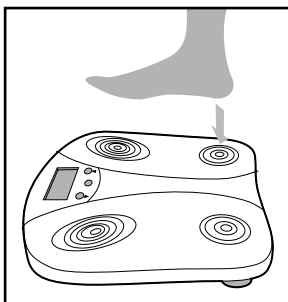
3. The gender indicator is then enabled, for male and for female. Press either or to change gender entry. Press SET to confirm.

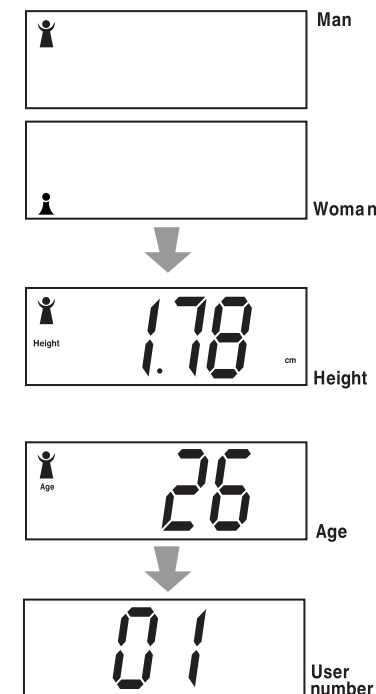
4. The unit is then ready to accept your height. Press or to change to your correct height. The range is from 100cm - 220cm (1m-2.2m) or 3'4"-7'2". Press SET to confirm.

5. The unit is ready to accept your age. Press or to increase or decrease the displayed number of your age. The range is 10 to 100. Press SET to confirm.

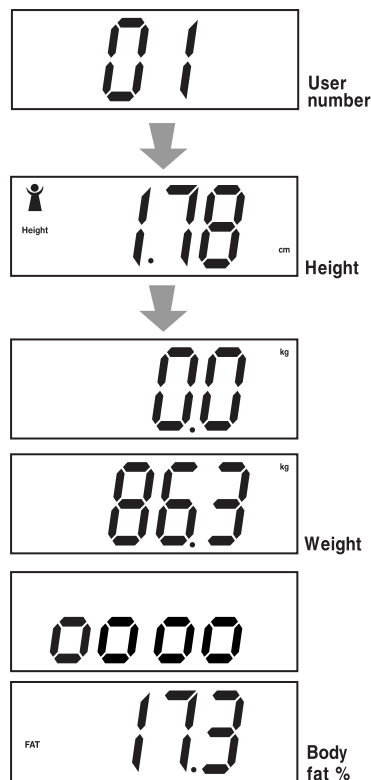
6. The user number will be displayed again. Press the SET key again until "0" is shown in the display. The unit is now ready for measuring your body fat.

**Note:**

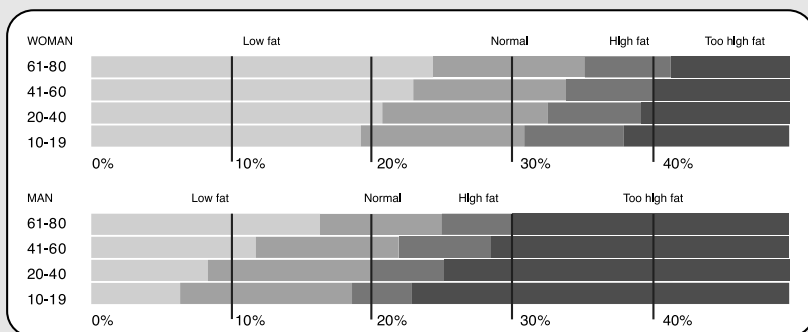
1. If there is a delay of more than about 15 seconds when entering personal data, the unit will automatically turn off.
2. The user must complete the steps 1 - 6 above to ensure a valid data entry.
3. If there is no pre-entered information, the unit assumes the data to be male, 1.80 m in height and age 30.

## ENGLISH

**MEASURING BODY FAT PERCENTAGE**

1. Step on the scale with bare feet and the display will show progressively the user numbers from 01 to 04. When the desired user number is displayed step off the scale immediately.
2. The programmed information about gender, height and age for the selected user will be displayed. The display will return to "0.0".
3. Stand on the unit again with bare feet placing your soles and heels on the rings and wait until your weight is displayed.
4. Continue standing on the scale and the display will then show "0000".
5. After a few seconds, the body fat percentage is displayed.

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